

# List of CDS Related Guidance


Processing Requirements for Thermolabile Flexible Endoscopes



# CDS Related Guidance

## Processing Requirements for Thermolabile Flexible Endoscopes

Denmark

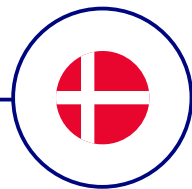

Procedure Room	Decontamination Area				
Precleaning	Manual Cleaning*	Manual Disinfection or EWD	Drying/Storage/ Sterilization	Remarks	Source
Immediately after finishing patient's procedure  By wiping the outer surface of the insertion tube and using <b>water</b> to flush/rinse all channels	Mainly enzymatic cleaners (but also mild alkaline cleaners)	EWD processing is state of the art  <b>PAA</b> in preference to GDA	Use of drying/storage cabinets	Flexible cystoscopes and bronchoscopes classified as critical medical devices require sterilization (if applicable, H <sub>2</sub> O <sub>2</sub> sterilization) after all manual (pre-) cleaning steps and EWD	Statens Serum Institut; Nationale Infektionshygiejniske Retningslinjerl. For Genbehandling af Fleksible Endoskoper (Version 6.2, <b>2021</b> ) <a href="https://hygiejne.ssi.dk/NIRendoskoper">https://hygiejne.ssi.dk/NIRendoskoper</a>

<b>Legend:</b> EWD – endoscope washer-disinfector (chemo-thermal mechanical/ automatic cleaning and disinfection) GDA – glutardialdehyde	PAA – peracetic acid ClO <sub>2</sub> – chlorine dioxide LCS – liquid chemical sterilization EtO – ethylene oxide	LTSF – low-temperature steam formaldehyde H <sub>2</sub> O <sub>2</sub> – hydrogen peroxide OR – operating room EHEF – European Hygiene Experts Forum
--	--	--

\*Mandatory, even in the case of EWD processing

# List of CDS Related Guidance – Denmark

As medical knowledge is constantly growing, technical modifications or changes of the product design, product specifications, accessories and service offerings may be required. This quick reference guide is a mere summary of some important procedure steps and/or product information. It is not meant to replace the instructions for use. Any user of this product must at all times observe all mandatory information for the product, found, in particular, on the labels and the instructions for use.



**OLYMPUS EUROPA SE & CO. KG**  
Postbox 10 49 08, 20034 Hamburg, Germany  
Wendenstrasse 20, 20097 Hamburg, Germany  
Phone: +49 40 23773-0, Fax: +49 40 233765  
[www.olympus-europa.com](http://www.olympus-europa.com)