


Single Use Injector

# NM-221C-0427

**High Kink-resistant Sheath and Retractable Needle for Safe and Smooth Injection**

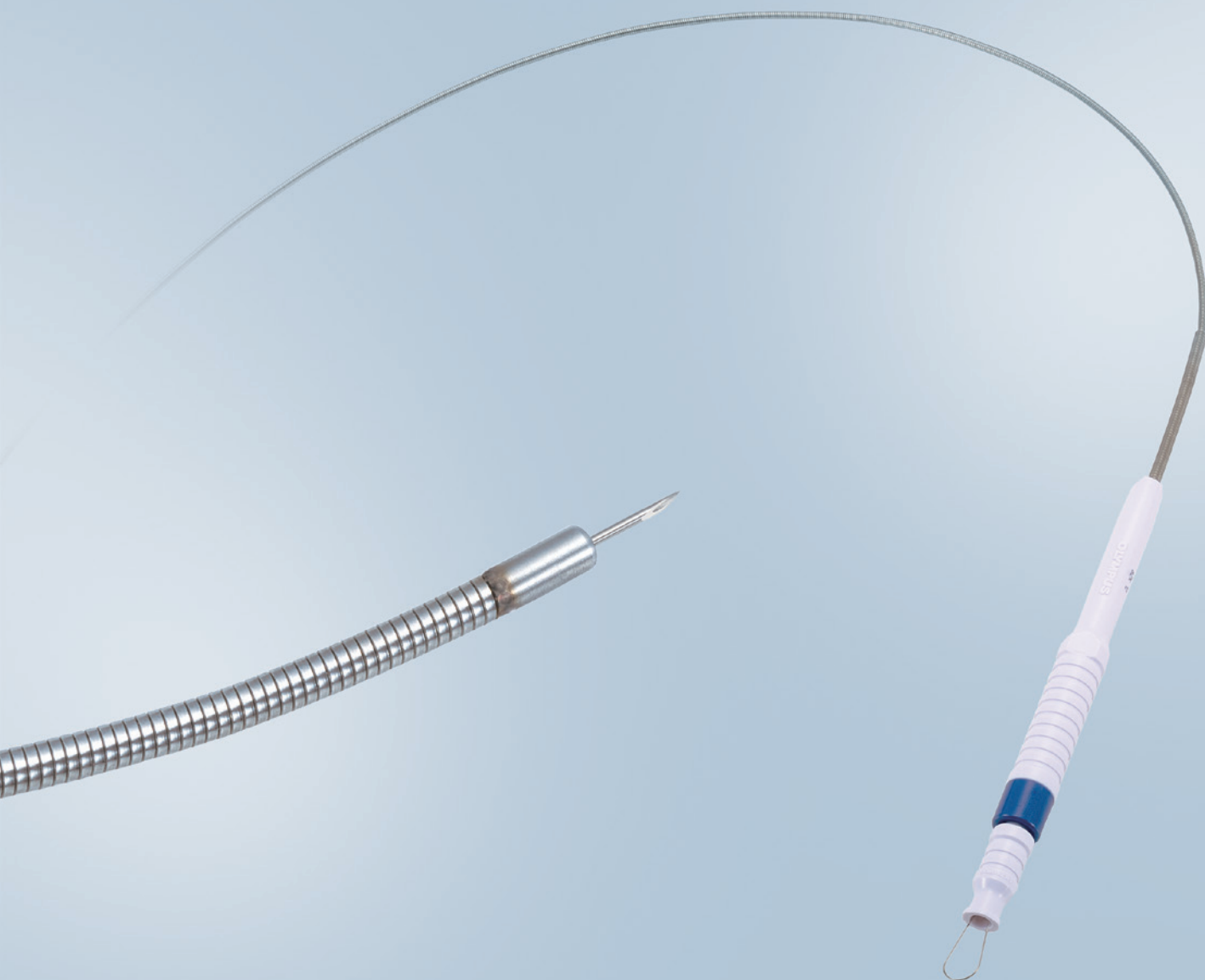

# Single Use Injector NM-221C-0427

## Reliable Design for Multiple Needle Injection into Bladder Wall

### Metal Coil Sheath Design

Metal coil sheath design is highly kink-resistant which allows for multiple injections of the needle into the bladder wall. The sharp, durable needle grind profile also supports smooth multiple injections.



### Retractable Needle with Acoustic Feedback

The retractable needle design protects the scope channel from needle puncture. In addition, there is an acoustic feedback so that the physician can easily verify that the needle has been retracted into the sheath. This also prevents mishandling and scope damage by unintended needle advancement.


### Ready-to-Use Sterilized Package

Single-use design eliminates the inconvenience of reprocessing, and helps mitigate the risk of infection. The ready-to-use sterilized package allows the physician to use the device immediately.

Specifications	
Model	NM-221C-0427
Maximum Insertion Portion Diameter	1.84 mm
Working Length	971 mm
Needle Diameter	27 G
Needle Length	4 mm
Unit	1 sheath & 1 needle / Box

Compatible Olympus Endoscopes	
Model and Length	CYF; working length up to 380 mm
Minimum Channel Diameter	2 mm

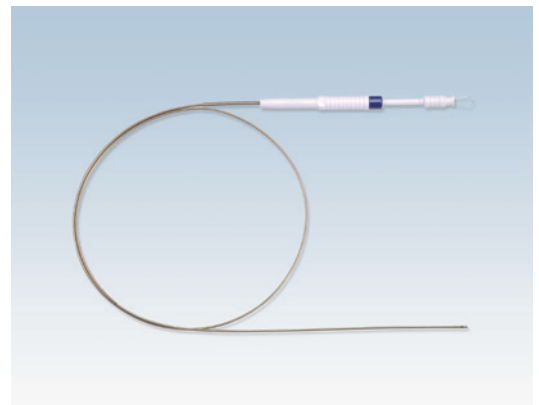

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

**OLYMPUS EUROPA SE & CO. KG**  
Postbox 10 49 08, 20034 Hamburg, Germany  
Wendenstrasse 14-18, 20097 Hamburg, Germany  
Phone: +49 40 23773-0, Fax: +49 40 233765  
www.olympus-europa.com