


BiCOAG™ Hemostasis Probe

BiCOAG™

Optimized tissue-to-electrode contact for effective and efficient hemostasis

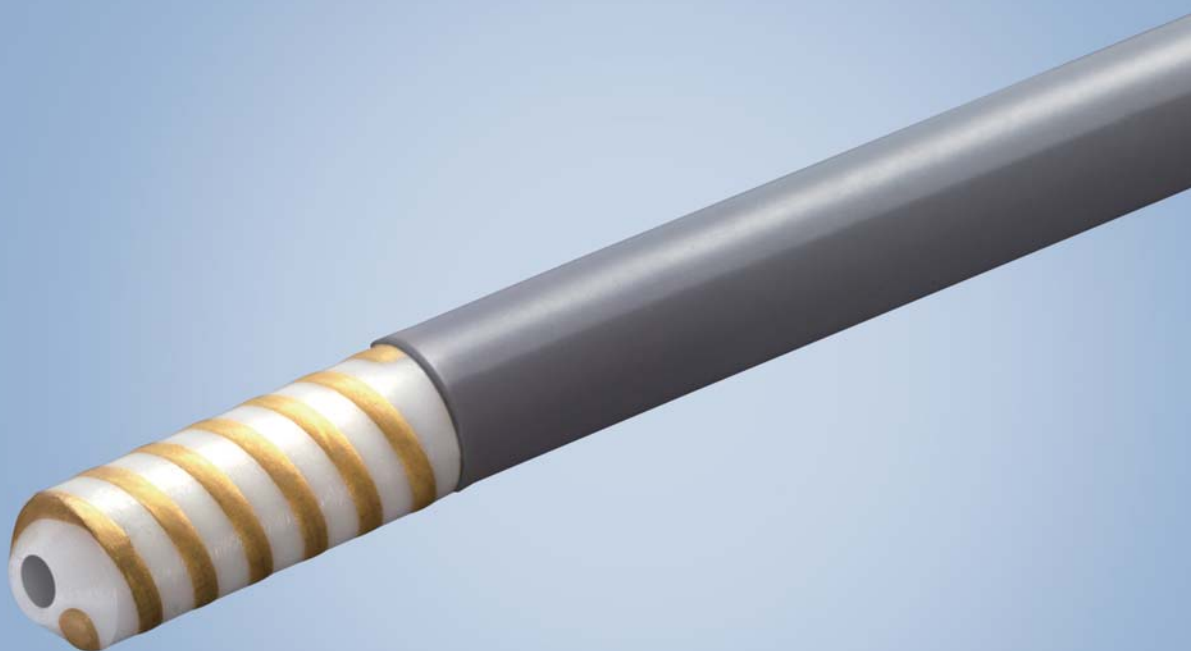


The single use BiCOAG™ hemostasis probe provides reliable hemostasis through bipolar coagulation technology.




Bipolar Spiral Tip

- Coagulation at virtually any angle

Probe Side View

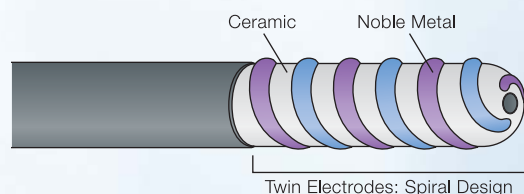

Probe End View

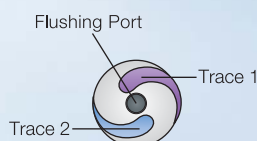

Optimal Tubing Stiffness

- Reduces kinking
- Transmits pressure to bleeding site
- Passes easily through the endoscope channel

Flushing Port


- Central lumen for irrigation

Choice of Connectors

- Coaxial Connector
- Fixed Pin Connector



Coaxial Connector

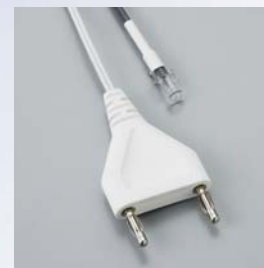

Fixed Pin Connector

Specifications

Model	Packed Quantity	Maximum Probe Tip Outer Diameter	Connector Type	Flushing Port	Working Length	Compatible Channel Inner Diameter
CD-B610LA	1 piece/box	10 Fr	Coaxial	Available	3500 mm	3.7 mm
CD-B612LA			Fixed Pin			
CD-B620LA		7 Fr	Coaxial			2.8 mm
CD-B622LA			Fixed Pin			

Electro Surgical Generator / Peristaltic Pump ESG-100 / AFU-100

Compact and easy-to-use, the ESG-100 features a variety of monopolar and bipolar modes for cutting and coagulation, and is designed to work together with the peristaltic pump AFU-100.



ESG-100



AFU-100

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.