


Multi-3 V Plus™

Single Use 3-Lumen Extraction Balloon V



Excellent balloon durability combined with efficient stone extraction capability provides choice and flexibility for ERCP physicians

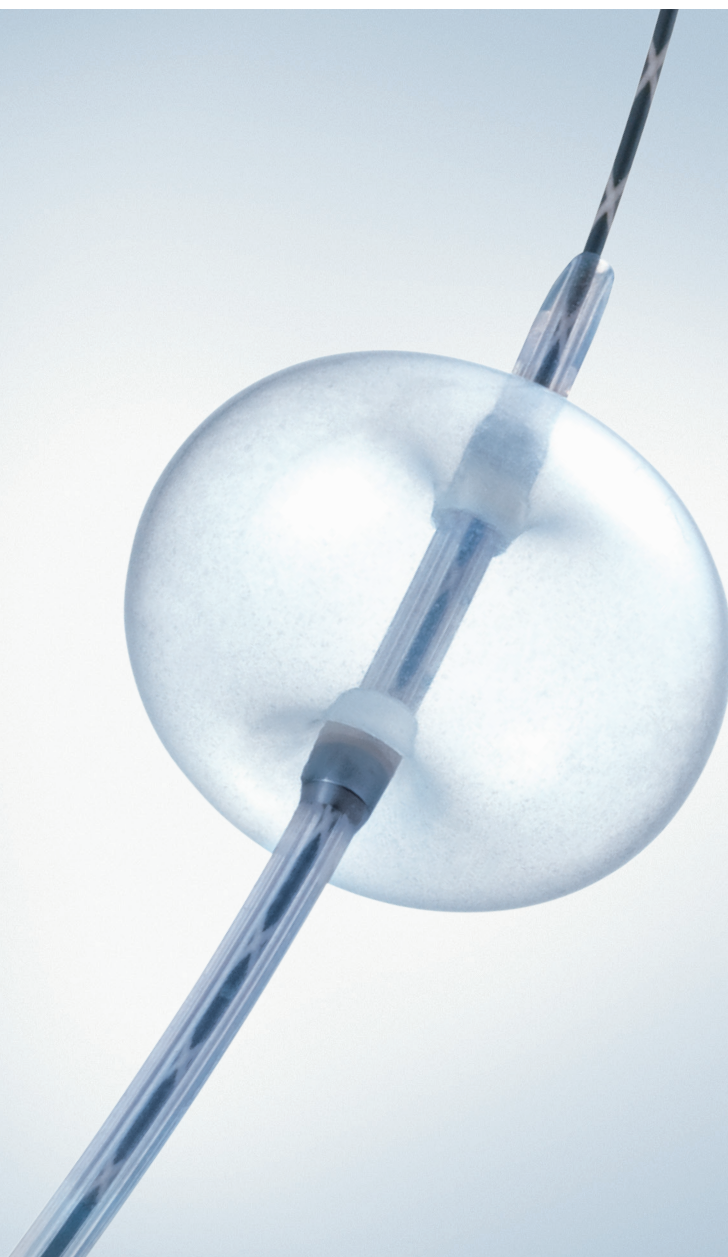


Excellent balloon durability supports repeated stone extraction.

Excellent balloon durability

Thick, high-quality latex increases the durability and reduces the possibility of the balloon bursting during stone extraction


Multiple balloon sizing

Balloons can be inflated to one of three diameters: 8.5 mm, 11.5 mm, and 15 mm or 15 mm, 18 mm, and 20 mm. Each balloon's size can be easily adjusted to suit the anatomical conditions of each case.

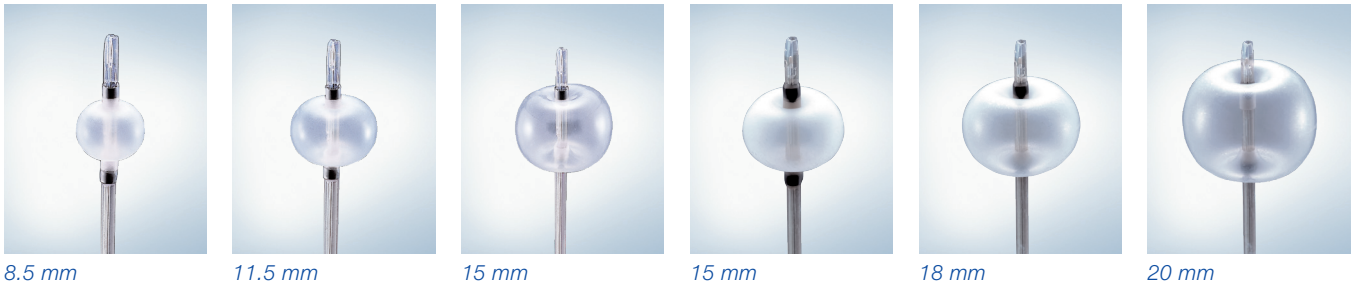


Choice of distal wireguided or over the wire design

Distal wireguided option facilitates extraction of multiple stones with the guidewire still in the bile duct.

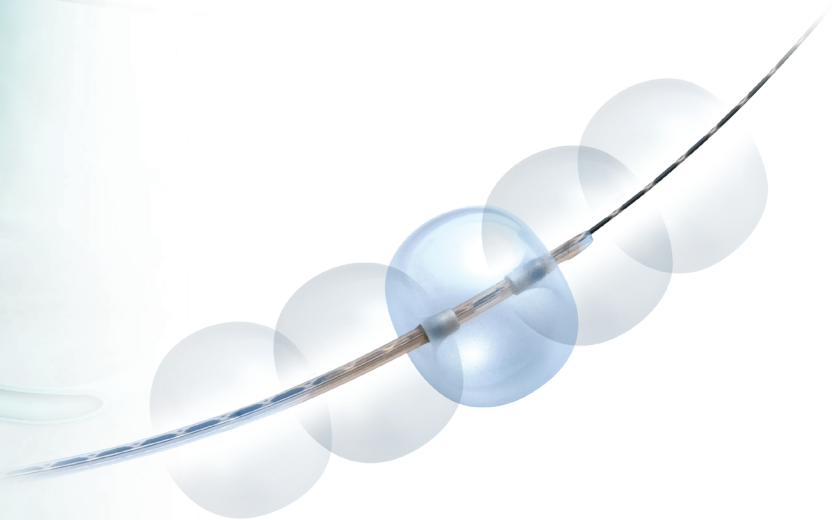


Three pre-measured syringes for reliable inflation

Clearly marked syringes enable precise inflation of each balloon quickly without having to visually check the markings on the syringe.

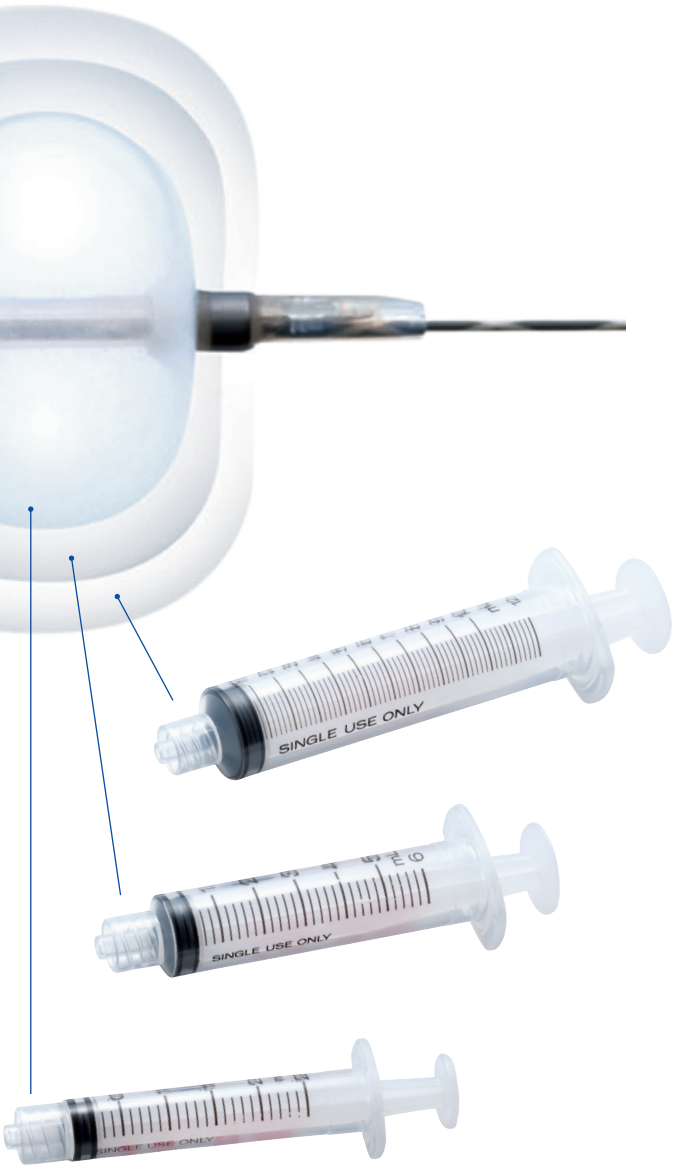

Specifications				
HF System Components				
Model	B-V232P-A/B	B-V242Q-A/B	B-V432P-A/B	B-V442Q-A/B
Maximum Balloon Diameter	15 mm	15 mm	20 mm	20 mm
Sheath Diameter		Distal end: 5.5 Fr	Proximal end: 7 Fr	
Sheath Design	Over the wire	Distal wireguided	Over the wire	Distal wireguided
Working Length	1900 mm	1950 mm	1900 mm	1950 mm
Minimum Channel Size	2.8 mm	3.2 mm	3.2 mm	3.2 mm
Compatible Guidewire	0.89 mm (0.035 inch)			
Radiopaque Band	One band at the distal end and one band at the proximal end of balloon			
Syringe	Three syringes are contained in the balloon package			
	8.5 mm/11.5 mm/15 mm	8.5 mm/11.5 mm/15 mm	15 mm/18 mm/20 mm	15 mm/18 mm/20 mm
Lumen	Triple Lumen			
Injection Port	A: Above type B: Below type			

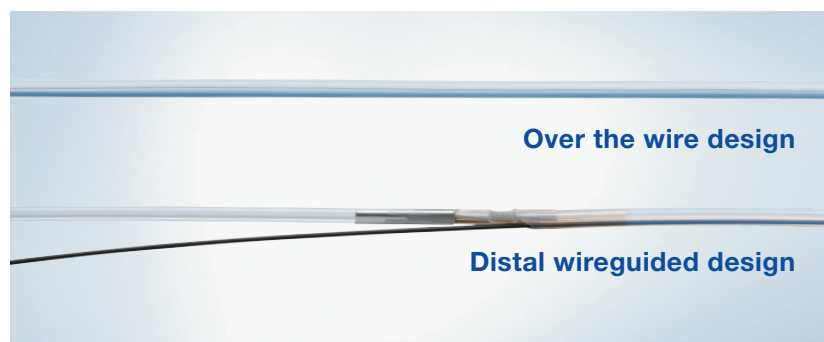

Olympus EndoTherapy

Dedicated Support beyond Products

At Olympus, we strive to be more than just a medical-equipment provider to our customers. We provide end-to-end support, from the purchasing process to procedure and reprocessing services, to build a relationship of trust.



ERCP Showroom Page

Enter the world of ERCP on our Olympus website for a full overview of our product portfolio, valuable video content, information on training and much more.

> www.olympus.eu/ercp-procedure



Training

Register here to improve procedure skills and the safe and effective use of Olympus products.

> www.olympus.eu/training



Portfolio

Browse through the rich portfolio of EndoTherapy devices here.

> www.olympus.eu/et-catalog



Service

Contact your Olympus sales representative for individual support now.

As medical knowledge is constantly growing, technical modifications or changes of the product design, product specifications, accessories and service offerings may be required.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG
Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 20, 20097 Hamburg, Germany
Phone: +49 40 23773-0, Fax: +49 40 233765
www.olympus-europa.com