

Efficiency and Choice


DIEGO ELITE – the Olympus Rhinology Shaver System



DIEGO ELITE Multidebrider

Efficiency – Aiming to Improve Daily Surgery

The DIEGO ELITE Multidebrider is a multipurpose ENT shaver system designed to offer surgeons improved choice on how to perform surgeries while at the same time offering features aimed at improving surgical efficiency.

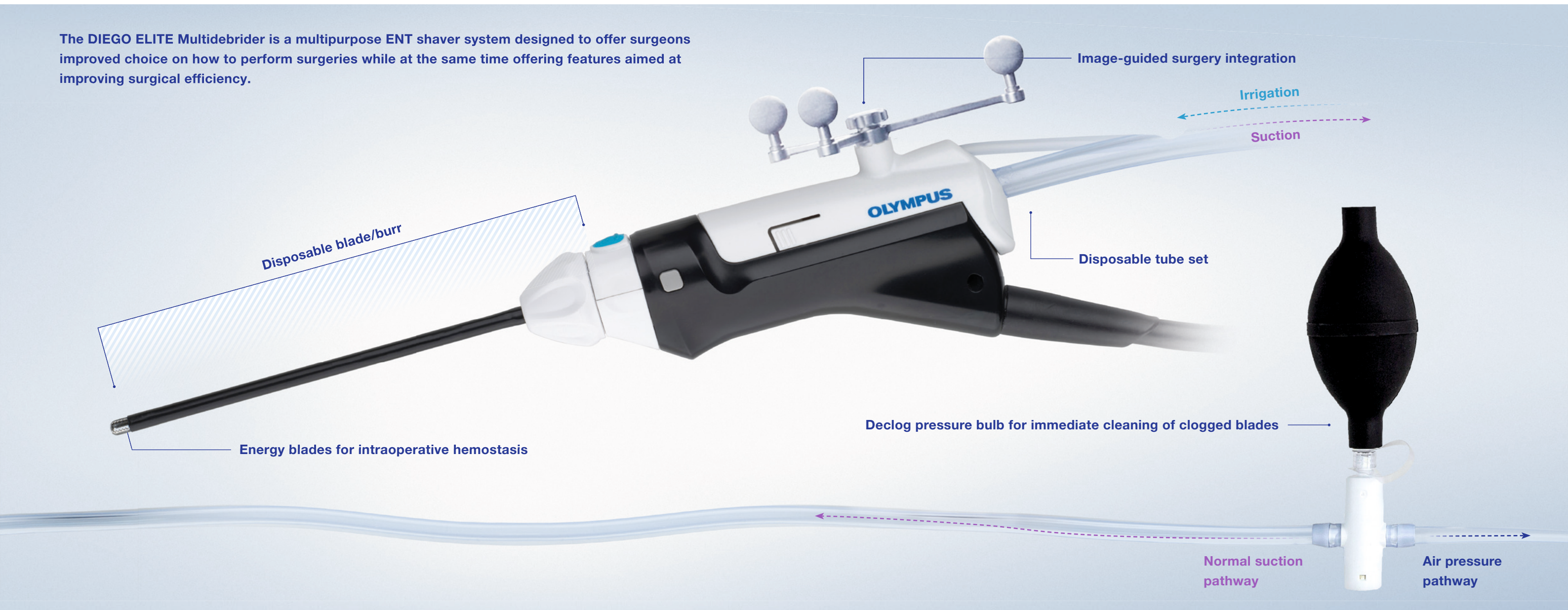

Full System Overview

Touchscreen Console

- Intuitive and simple handling.
- Blade ID recognition with predefined optimal setting.
- Dedicated suction unit.



Reusable Handpieces

- Straight and angled option.
- Comfortable and ergonomic handling.
- Open design for safer cleaning. Autoclavable.



Suction Tubing

- Completely disposable suction pathway.
- Declog suction tubing for immediate cleaning of clogged blades.



Toggle Stop

- Improved efficiency and safety. One touch to have the blade's window fully open/closed, avoiding undesired tissue engagement or having a suction in your hand.





Choice — Multiple Ways to Perform Surgery

Blade Portfolio

The DIEGO ELITE Multidebrider offers the surgeon multiple ways to perform surgery, including a choice of two handle types, comprehensive portfolios of blades and burrs as well as the possibility to choose between standard or declog tube set.


Standard Blades

Standard soft tissue removal in the nasal anatomy.



Energy Blades

Bipolar and monopolar energy blades for immediate hemostasis during surgery.

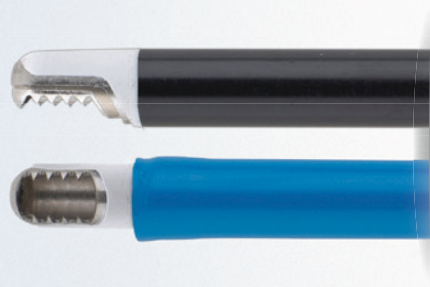

Rotatable Tip Blades

Allows the surgeon to reach different areas of the nasal anatomy without twisting the handpiece.



Type A Blades

For advanced resection.




Angled Blades

A variety of angles including 0°, 15°, 60°, 75° and 90°.



Turbinate Blades

2 mm turbinate blades for submucosal resection also available with bipolar option.

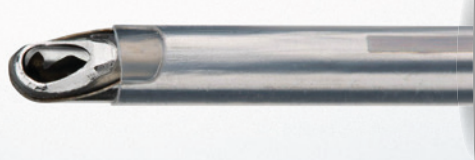

Beyond Functional Endoscopic Sinus Surgery (FESS)

Tonsil and Adenoid Blades

A wide range of blades are available for tonsil and adenoid resection. Malleable, bipolar and monopolar blades for cauterization and stabilization of the adenoid bed.



Burr Portfolio

The DIEGO ELITE system offers a comprehensive selection of burrs for advanced FESS and skull base procedures

Diamond Bullet Burr

Diamond bullet burr portfolio for gentle removal of bone and cartilage structures. Available in different angulations.



Diamond Ball Burr

Diamond removal of delicate bone and cartilage structures. Available in different angulations.



Angled Burrs

Full burr portfolio available in a variety of angles including 15°, 60° and 70°.



Diamond Taper Burr

Diamond taper burr portfolio for enlarging openings and for bone and cartilage removal, especially designed to work from back to front increasing visualization. Available in different angulations.



Fluted Barrel Burr

Cutting burr portfolio for advanced removal of bone and cartilage structures. Available in different angulations.



Laryngeal Blades

For the removal of laryngeal lesions and tumor debulking. Available as 2.9 mm and 4 mm with a length of 22 cm for subglottic / upper tracheal area or 27 cm for glottic / supraglottic applications.



Did You Know ...

Olympus Offers a Full See-and-Treat Portfolio for Office- and OR-Based Rhinology Procedures

Imaging

Treatment



As medical knowledge is constantly growing, technical modifications or changes of the product design, product specifications, accessories and service offerings may be required.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG
Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 20, 20097 Hamburg, Germany
Phone: +49 40 23773-0, Fax: +49 40 233765
www.olympus-europa.com